

All data taken at the Pacific Northwest National Laboratory

FTS Operators: Tracy J. Baker

Data Analysis: Tracy J. Baker

Composite spectrum for: 1-Octanol

- First Column: Position in wavenumber (cm^{-1})
- Second column: Real refractive index $n(\tilde{\nu})$ (dispersion index)
- Third column: Imaginary refractive index, $k(\tilde{\nu})$ (absorption index per unit length in centimeters)

Where the complex refractive index $\hat{n} = n(\tilde{\nu}) + ik(\tilde{\nu})$

Following Bertie (in the references below) we define the absorbance as $A = -\log_{10}(I/I_0)$ and the linear absorption coefficient $K = A/d$, where d is the path length. The connection between the imaginary refractive index and the absorbance coefficient arises from the following: $2.303K = 4\pi\tilde{\nu}k$

See the following references for a detailed description of terms and units:

- 1) Bertie, J. E., Zhang, S. L., Eysel, H. H., Baluja, S., & Ahmed, M. K. (1993). Infrared Intensities of Liquids XI: Infrared Refractive Indices from 8000 to 2 cm^{-1} , Absolute Integrated Intensities, and Dipole Moment Derivatives of Methanol at 25°C . *Appl. Spec.*, 47(8), 1100-1114 doi:10.1366/0003702934067973
- 2) Bertie, J. E., Zhang, S. L., & Keefe, C. D. (1995). Measurement and use of absolute infrared absorption intensities of neat liquids. *Vibrational Spectroscopy*, 8(2), 215-229. doi:10.1016/0924-2031(94)00038-i

Sample:

- Chemical name, formula and CAS number: 1-Octanol, $\text{C}_8\text{H}_{18}\text{O}$, [111-87-5]
- IUPAC name: Octan-1-ol
- Synonyms: Alcohol C_8 ; Capryl alcohol; Octyl alcohol
- Physical properties: FW = 130.23 g/mole; mp = -15°C ; bp = $194\text{--}195^\circ\text{C}$; $\rho = 0.827\text{ g/cm}^3$
- Supplier and stated purity: TCI, 99% (Lot # OBWRM-PI)
- Temperature of sample: 25°C ($\pm 1^\circ\text{C}$)
- Individual samples were measured at the following path lengths: MIR: 4.15, 5.81, 9.08, 13.7, 17.2, 35.0, 46.6, 86.8, 216 and 508 micrometers (μm); NIR: 86.6, 216, 510, 920 and 4140 μm . Final data are a composite of these spectra.
- Sample cell window material: MIR = potassium bromide (KBr) except potassium chloride (KCl) for the 216 and 508 μm cells; NIR = KCl except KBr for the 86.6 and 920 μm cells.
- Preparation: None

NIR Instrument Parameters:

- Bruker Vertex 70, purged with UHP nitrogen
- Spectral range: 10,000 to $3,000\text{ cm}^{-1}$ (1.0 to 3.33 microns)
- NIR source: Quartz tungsten bulb
- Beamsplitter: Broadband potassium bromide (KBr)
- Detector: DLTGS at room temperature
- Aperture: 3 mm
- Folding limits: 31597.6 to 0 cm^{-1}


MIR Instrument Parameters:

- Tensor II with sample stage rotated 6 degrees, purged with UHP nitrogen
- Spectral range: 7,800 to 400 cm^{-1} (1.282 to 25 microns)
- IR source: Silicon carbide glow bar
- Beamsplitter: Potassium bromide (KBr)
- Detector: DTGS at room temperature
- Aperture: 3 mm
- Folding limits: 11664.15 to 0 cm^{-1}

NIR/MIR Instrument Parameters:

- Instrument resolution: 2.0 cm^{-1}
- Number of interferograms averaged per single channel spectrum: 128
- Apodization: Norton-Beer, Medium
- Phase correction: Mertz
- Scanner velocity: 10 kHz; 7.5 kHz (MIR)
- Interferogram zerofill: 4x
- Spectral interval after zerofilling: 0.4822 cm^{-1}

a)

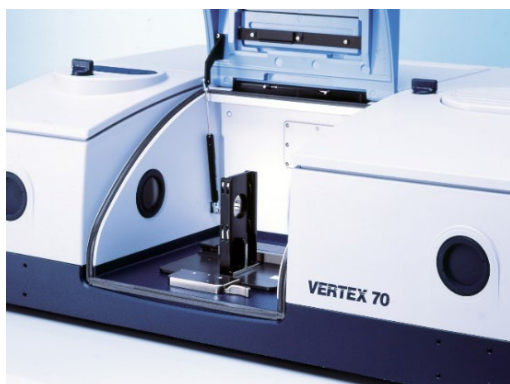

b)



Figure 1: The Bruker Vertex 70 FTIR (a) and Tensor II (b).

Measured Refractive Index:

The refractive index for 1-Octanol was measured at $25\text{ }^{\circ}\text{C}$ using an Atago model DR-M2/1550 Abbe refractometer. Notch filters were employed in front of a white light source to make measurements at multiple wavelengths. An infrared viewer from Atago was used to detect signal at 1550 nm . The temperature was controlled to match that in the sample compartment of the FTIR using a heated circulating bath.

480 nm:	$n = 1.4323$	486 nm:	$n = 1.4319$	546 nm:	$n = 1.4286$
589 nm:	$n = 1.4268$	644 nm:	$n = 1.425$	656 nm:	$n = 1.4246$
1550 nm:	$n = 1.4097$				

The refractive index, n , vs. wavelength in microns, λ , was fit to an equation similar to that of Sellmeier:

$$n(\lambda) = \{a + b/(\lambda^2 - c)\}^{1/2}$$

The resulting best-fit equation was used to find the refractive index at the highest energy data points in our experimental spectra. For 1-Octanol, the results were

$$\begin{aligned} n(7,800\text{ cm}^{-1}) &= 1.4118 \text{ at } 25\text{ }^{\circ}\text{C} \text{ for MIR data and} \\ n(10,000\text{ cm}^{-1}) &= 1.4155 \text{ at } 25\text{ }^{\circ}\text{C} \text{ for NIR and merged data.} \end{aligned}$$

Post Processing and Related Parameters:

For the MIR, a composite spectrum was created from 10 absorbance spectra (base-10) taken at 10 path lengths: 4.15, 5.81, 9.08, 13.7, 17.2, 35.0, 46.6, 86.8, 216 and 508 micrometers (μm). These data were collected with the sample stage rotated by -6° to minimize artifacts resulting from back-reflection into spectrometer (see Johnson et al., *Appl. Spectrosc.*, 76(5) 620-624, 2021). For the NIR, a composite spectrum was created from 5 absorbance spectra (base-10) taken at 5 path lengths: 86.6, 216, 510, 920 and 4140 μm . At each path length several spectra were measured and the results averaged for better signal to noise. The measured cell lengths were adjusted using Beer's law plots in which the NIR and MIR data were analyzed independently.

- 1) The imaginary part of the refractive index, or k vector, was determined for each absorbance file as per Bertie's program "RNJ46A" (see reference above). This takes into account the reflective losses due to the KBr and/or KCl windows.
- 2) A composite k vector is created via a classical, weighted, linear, least squares fit using the output files of program "RNJ46A": Intercept=0, slope is fitted, individual absorbance values weighted by T^2 (transmission squared), all absorbance values ≥ 2.5 are given zero weight. For the MIR, eleven composite vectors were created and merged by hand.
 - a) The first k vector used the results from the 508 μm cell. This k vector determined the final values for the range from 7800 to 3560 cm^{-1} .
 - b) The second k vector used the results from the 9.08 through 35.0 μm cells. This k vector determined the final values for the range from 3560 to 3450 cm^{-1} and 2845 to 2600 cm^{-1} .
 - c) The third k vector used the results from the 5.81, 13.7 and 17.2 μm cells. This k vector determined the final values for the range from 3450 to 3250 cm^{-1} .
 - d) The fourth k vector used the results from the 5.81 and 13.7 through 35.0 μm cells. This k vector determined the final values for the range from 3250 to 2970 cm^{-1} .
 - e) The fifth k vector used the results from the 5.81 and 9.08 μm cells. This k vector determined the final values for the range from 2970 to 2940 cm^{-1} and 2915 to 2848 cm^{-1} .
 - f) The sixth k vector used the results from the 5.81 μm cell. This k vector determined the final values for the range from 2940 to 2915 cm^{-1} .
 - g) The seventh k vector used the results from the 216 and 508 μm cells. This k vector determined the final values for the range from 2600 to 1800 cm^{-1} .
 - h) The eighth k vector used the results from the 216 μm cell. This k vector determined the final values for the range from 1800 to 1500 cm^{-1} .
 - i) The ninth k vector used the results from the 13.7 through 46.6 μm cells. This k vector determined the final values for the range from 1500 to 1020 cm^{-1} .
 - j) The tenth k vector used the results from the 17.2 through 86.8 μm cells. This k vector determined the final values for the range from 1020 to 720 cm^{-1} .
 - k) The eleventh k vector used the results from the 35.0 through 86.8 μm cells. This k vector determined the final values for the range from 720 to 370 cm^{-1} .
- 3) A frequency correction was applied to the resulting composite MIR k vector.
 - a) Frequency correction (already applied): $\tilde{\nu}(\text{corrected}) = [\tilde{\nu}(\text{instrument}) * 0.99988 + 0.025]$ as determined by comparing measured atmospheric spectral lines (H_2O and CO_2) to values from the Northwest Infrared Spectral Library Database.
- 4) For the NIR, four composite vectors were created and merged by hand.
 - a) The first k vector used the results from the 4140 μm cell. This k vector determined the final values for the range from 10,000 to 4500 cm^{-1} .
 - b) The second k vector used the results from the 216 and 510 μm cells. This k vector determined the final values for the range from 4500 to 4250 cm^{-1} and 3650 to 3600 cm^{-1} .
 - c) The third k vector used the results from the 510 and 920 μm cells. This k vector determined the final values for the range from 4250 to 3650 cm^{-1} .
 - d) The fourth k vector used the results from the 86.6 through 510 μm cells. This k vector determined the final values for the range from 3600 to 370 cm^{-1} .
- 5) The resulting composite NIR k vector and the refractive index at 10,000 cm^{-1} were used to create the real or n vector using the Kramers-Kronig relation, as per Bertie's program "LZZKTB."
 - a) Frequency correction (already applied): $\tilde{\nu}(\text{corrected}) = [\tilde{\nu}(\text{instrument}) * 0.99998 - 0.0005]$ as

determined by comparing measured atmospheric spectral lines (H_2O and CO_2) to values from the Northwest Infrared Spectral Library Database.

- 6) The MIR data were mapped onto the NIR x-axis using an interpolation routine, i.e. the Make Compatible command in OPUS 5.5. Then the composite MIR and NIR k vectors were merged to generate a final composite k vector across the entire spectral range. The NIR data were used exclusively above 3718 cm^{-1} , and only the MIR data were used below 3574 cm^{-1} . A weighted average, with the weight of the MIR vector increasing linearly from 0 to 100% between 3718 and 3574 cm^{-1} was used in the overlapping spectral region. The resulting composite k vector and the refractive index at $10,000\text{ cm}^{-1}$ were used to create the final n vector using the Kramers-Kronig relation, as per Bertie's program "LZZKTB."

Photograph of Sample 1-Octanol:

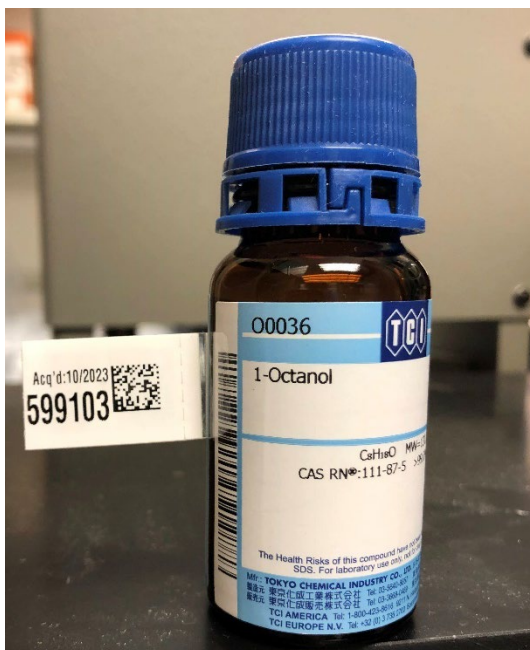


Figure 2: 1-Octanol in TCI container for NIR and MIR measurements.