

All data taken at the Pacific Northwest National Laboratory
FTS Operators: Gerry Gonzalez, Ashley M. Oeck and John S. Loring
Data Analysis: Russell G. Tonkyn

Composite spectrum for: Di-n-octyl phthalate

- First Column: Position in wavenumber (cm^{-1})
- Second column: Real refractive index $n(\tilde{\nu})$ (dispersion index)
- Third column: Imaginary refractive index, $k(\tilde{\nu})$ (absorption index per unit length in centimeters)

Where the complex refractive index $\hat{n} = n(\tilde{\nu}) + ik(\tilde{\nu})$

Following Bertie (in the references below) we define the absorbance as $A = -\log_{10}(I/I_0)$ and the linear absorption coefficient $K = A/d$, where d is the path length. The connection between the imaginary refractive index and the absorbance coefficient arises from the following: $2.303K = 4\pi\tilde{\nu}k$

See the following references for a detailed description of terms and units:

- 1) Bertie, J. E., Zhang, S. L., Eysel, H. H., Baluja, S., & Ahmed, M. K. (1993). Infrared Intensities of Liquids XI: Infrared Refractive Indices from 8000 to 2 cm^{-1} , Absolute Integrated Intensities, and Dipole Moment Derivatives of Methanol at 25°C . *Applied Spectroscopy*, 47(8), 1100-1114. doi:10.1366/0003702934067973
- 2) Bertie, J. E., Zhang, S. L., & Keefe, C. D. (1995). Measurement and use of absolute infrared absorption intensities of neat liquids. *Vibrational Spectroscopy*, 8(2), 215-229. doi:10.1016/0924-2031(94)00038-i

Sample:

- Chemical name, formula and CAS number: Di-n-octyl phthalate, $\text{C}_{24}\text{H}_{38}\text{O}_4$, [117-84-0]
- IUPAC name: Dioctyl benzene-1,2-dicarboxylate
- Synonyms: Dioctyl phthalate; Phthalic acid di-n-octyl ester; DNOP; n-Octyl phthalate
- Physical properties: FW = 390.56 g/mole; mp = -25°C ; bp 220°C ; $\rho = 0.980\text{ g/cm}^3$
- Supplier and stated purity: Aldrich, $\geq 98\%$ (Lot # BCBR9722V)
- Temperature of sample: 26°C ($\pm 1^\circ\text{C}$)
- Individual samples were measured at the following path lengths: MIR: 8.17, 10.7, 13.6, 38.2, 58.1, 102, 199 and $507\text{ }\mu\text{m}$; NIR: 123, 222, 504, 1225 and $3712\text{ }\mu\text{m}$. Final data are a composite of these spectra.
- Sample cell window material is potassium bromide (KBr) except the 504 (NIR) and 1225 (NIR) cell windows were potassium chloride (KCl).
- Preparation: None.

NIR Instrument Parameters:

- Bruker Vertex 70, purged with UHP nitrogen
- Spectral range: $10,000$ to $3,000\text{ cm}^{-1}$ (1.0 to 3.33 microns)
- NIR source: Quartz tungsten bulb
- Beamsplitter: Broadband Potassium bromide (KBr)
- Detector: DLTGS at room temperature
- Aperture: 3 mm
- Folding limits: 15801 to 0 cm^{-1}


MIR Instrument Parameters:

- Tensor 37, purged with UHP nitrogen
- Spectral range: 7500 to 400 cm^{-1} (1.33 to 25 microns)
- IR source: Silicon carbide glow bar
- Beamsplitter: Broadband Potassium bromide (KBr)
- Detector: DLTGS at room temperature
- Aperture: 3 mm
- Folding limits: 15802 to 0 cm^{-1}

NIR/MIR Instrument Parameters:

- Instrument resolution: 2.0 cm^{-1}
- Number of interferograms averaged per single channel spectrum: 128
- Apodization: Norton-Beer, Medium
- Phase correction: Mertz
- Scanner velocity: 10 kHz
- Folding limits: 15802 to 0 cm^{-1}
- Interferogram zerofill: 4x
- Spectral interval after zerofilling: 0.4823 cm^{-1}

a)

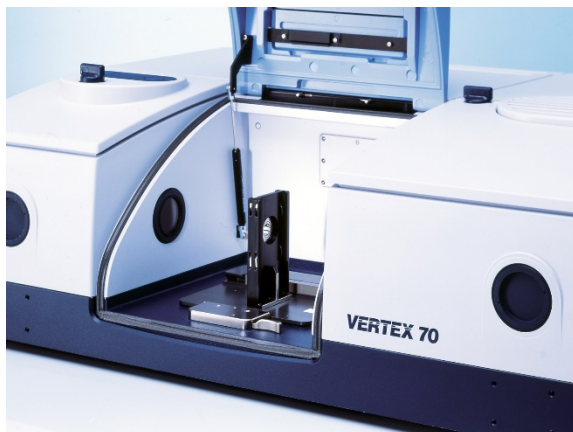

b)

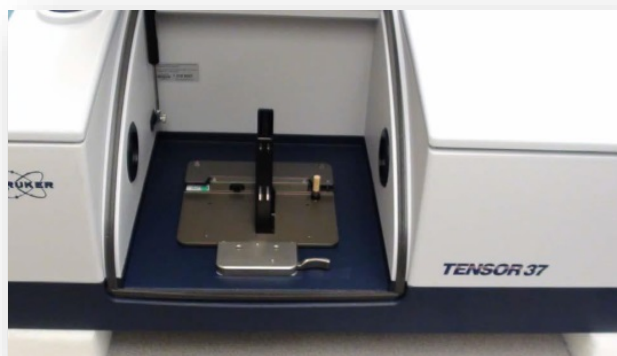


Figure 1: The Bruker Vertex 70 FTIR (a) and Tensor 37 FTIR (b).

Measured Refractive Index:

The refractive index for Di-n-octyl phthalate was measured at $25\text{ }^{\circ}\text{C}$ using an Atago model DR-M2/1550 Abbe refractometer. Notch filters were employed in front of a white light source to make measurements at multiple wavelengths. An infrared viewer from Atago was used to detect signal at 1550 nm . The temperature was controlled to match that in the sample compartment of the FTIR using a heated circulating bath.

| | | | | | |
|----------|--------------|---------|--------------|---------|--------------|
| 480 nm: | $n = 1.4905$ | 486 nm: | $n = 1.4900$ | 546 nm: | $n = 1.4846$ |
| 589 nm: | $n = 1.4817$ | 644 nm: | $n = 1.4791$ | 656 nm: | $n = 1.4783$ |
| 1550 nm: | $n = 1.4607$ | | | | |

The refractive index, n , vs. wavelength in microns, λ , was fit to an equation similar to that of Sellmeier:

$$n(\lambda) = \{a + b/(\lambda^2 - c)\}^{1/2}$$

The resulting best-fit equation was used to find the refractive index at the highest energy data point in our experimental spectra. For Di-n-octyl phthalate, the result was

$$\begin{aligned} n(7800\text{ cm}^{-1}) &= 1.4629 \text{ at } 25\text{ }^{\circ}\text{C} \text{ for MIR data and} \\ n(10,000\text{ cm}^{-1}) &= 1.4669 \text{ at } 25\text{ }^{\circ}\text{C} \text{ for NIR and merged data.} \end{aligned}$$

Post Processing and Related Parameters:

For the MIR, a composite spectrum was created from 8 absorbance spectra (base-10) taken at 8 path lengths: 8.17, 10.7, 13.6, 38.2, 58.1, 102, 199 and 507 micrometers (μm). For the NIR, a composite spectrum was created from 5 absorbance spectra (base-10) taken at 5 path lengths: 123, 222, 504, 1225 and 3712 μm . The same cells for the ~ 100 , 200 and 500 μm path lengths were used for both spectral ranges. At each path length several spectra were measured and the results averaged for better signal to noise. The measured cell lengths were adjusted using Beer's law plots.

- 1) The imaginary part of the refractive index, or k vector, was determined for each absorbance file as per Bertie's program "RNJ46A" (see reference above). This takes into account the reflective losses due to the KBr windows and the measured refractive index at 7800 cm^{-1} for the MIR data and $10,000\text{ cm}^{-1}$ for the NIR data.
- 2) A composite k vector is created via a classical, weighted, linear, least squares fit using the output files of program "RNJ46A": Intercept=0, slope is fitted, individual absorbance values weighted by T^2 (transmission squared), all absorbance values ≥ 2.5 are given zero weight. For the MIR, six composite vectors were created and merged by hand.
 - a) The first k vector used the results from the 200 and 507 μm cells. This k vector determined the final values for the range from 7800 to 4900 cm^{-1} .
 - b) The second k vector used the results from the 58 through 507 μm cells. This k vector determined the final values for the range from 4900 to 3300 cm^{-1} .
 - c) The third k vector used the results from the 8 through 507 μm cells. This k vector determined the final values for the range from 3300 to 2500 cm^{-1} .
 - d) The fourth k vector used the results from the 14 through 507 μm cells. This k vector determined the final values for the range from 2500 to 1900 cm^{-1} .
 - e) The fifth k vector used the results from the 8 through 58 μm cells. This k vector determined the final values for the range from 1900 to 650 cm^{-1} .
 - f) The sixth k vector used the results from the 58 and 507 μm cells. This k vector determined the final values for the range from 650 to 370 cm^{-1} .
- 3) A frequency correction was applied to the resulting composite MIR k vector.
 - a) Frequency correction (already applied): $\tilde{\nu}(\text{corrected}) = [\tilde{\nu}(\text{instrument}) * .99975 + .005115]$ as determined by comparing measured atmospheric spectral lines (H_2O and CO_2) to values from the Northwest Infrared Spectral Library Database.
- 4) For the NIR, two composite vectors were created and merged by hand.
 - a) The first k vector used the results from the 1225 and 3712 μm cells. This k vector determined the final values for the range from 10000 to 7500 cm^{-1} .
 - b) The second k vector used the results from the 123 through 3712 μm cells. This k vector determined the final values for the range from 7500 to 400 cm^{-1} .
- 5) A frequency correction was applied to the resulting composite NIR k vector.
 - a) Frequency correction (already applied): $\tilde{\nu}(\text{corrected}) = [\tilde{\nu}(\text{instrument}) * .999869 + .0158818]$ as determined by comparing measured atmospheric spectral lines (H_2O and CO_2) to values from the Northwest Infrared Spectral Library Database.
- 6) Finally, the MIR data were mapped onto the NIR x-axis using an interpolation routine, i.e. the Make Compatible command in OPUS 5.5. Then the composite MIR and NIR k vectors were merged to generate a final composite k vector across the entire spectral range. The NIR data were used exclusively above 5400 cm^{-1} , and only the MIR data were used below 4800 cm^{-1} . A weighted average, with the weight of the MIR vector increasing linearly from 0 to 100% between 5400 and 4800 cm^{-1} was used in the overlapping spectral region. The resulting composite k vector and the refractive index at $10,000\text{ cm}^{-1}$ were used to create the final n vector using the Kramers-Kronig relation, as per Bertie's program "LZZKTB."

Photograph of Sample Di-n-octyl phthalate:



Figure 2: Di-n-octyl phthalate in Aldrich container.